



Diffuser Attachment

For Mighty Scope ClearVue


Item # 26700-220-DIFF

Designed to be used with Aven's Mighty Scope ClearVue digital microscope. When in use, this attachment disperses the LED light from the Mighty Scope ClearVue for soft, even illumination!

Key Features

- Custom made attachment for the Mighty Scope ClearVue Digital Microscope
- Evenly disperses LED light from the Mighty Scope ClearVue for soft & even illumination
- Eliminates glares and hot spots caused by bright lights
- Ideal when examining shiny or reflective parts
- Easy to install and remove

**With
Diffuser**

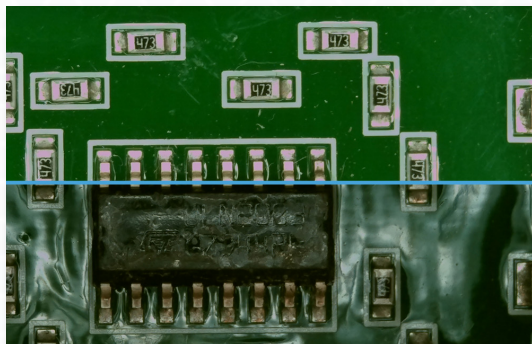

**Without
Diffuser**

**With
Diffuser**




**Without
Diffuser**

**Soft & Even
Illumination**

Eliminates glare
& hot spots from
shiny objects

**Easy to
Install**

Attaches to the
ClearVue with no
tools required

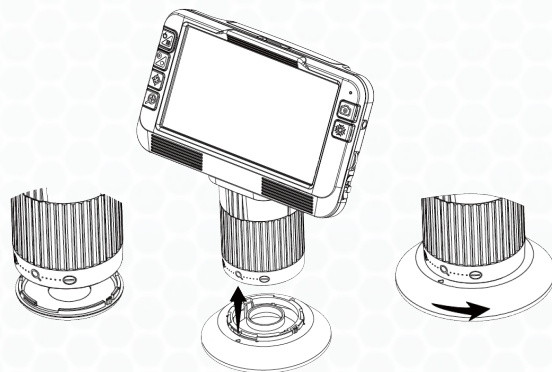

Easy Attachment

1. Remove the cap from the Mighty Scope ClearVue
2. Connect the diffuser attachment to the Mighty Scope ClearVue
3. Twist the diffuser attachment to the right to secure it

Diffuser Attachment For Mighty Scope ClearVue

Item #	26700-220-DIFF
Compatibility	Mighty Scope ClearVue Digital Microscopes
Material	Plastic
Recommended Working Distance	2.5 in - 3.5 in (when in use) (6.35 cm - 8.89 cm)
Dimensions	4 in x 4.75 in x 1.25 in 10.16 cm x 12.07 cm x 3.18 cm
Weight	0.05 lb